


Volk New Product Announcement


6 Mirror (G6) Gonio Lens



We are pleased to introduce the new 6 mirror lens - a truly great advancement for Gonioscopy. The 4 mirror lens has generally been the most popular lens for examining the anterior segment. The 6 mirror design provides comprehensive sweeping views of the anterior segment. No mirror gaps, no need to turn the lens. In essence a real 360° Gonioscopy examination.

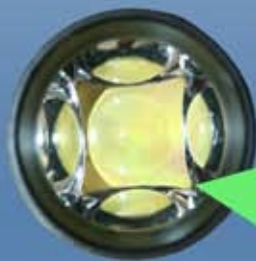

Note the mirror gaps in the 4 mirror; the lens needs to be turned to see the angle completely. Compare to the continuous image of the 6 mirror. This sweeping image will give a true 360° view of the angle.

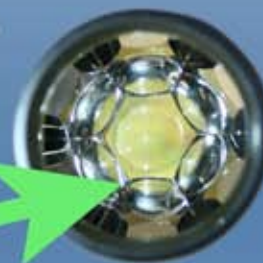

The method of examination remains the same. Two mirrors are aligned in the superior and inferior positions. The remaining four mirrors provide the continuously sweeping view through the nasal and temporal regions. In fact this method is a marked improvement. With a 1, 2, 3 or 4 mirror, the examiner had to note where he was geographically in the chamber angle and mirror location to move the lens or switch to another mirror. The 6 mirror with no gaps eliminates this need - one continuous view.

Lens	Type	Mirror Angles	Image Mag.	Laser Spot	Contact Diameter
VG6LNF	Ring	6 x 63°	1.0x	1.0x	8.4mm

See the new buck slip (ML-1520) for additional information on the features and benefits of this lens. Contact customer service for further information, pricing and availability.



7893 Enterprise Drive, Mentor, OH 44060.

• Tel: 800.345.8655 / 440.942.6161 • Fax: 440.942.2257 • Web: www.volk.com